EN

Manual

OptiSmile Water Flosser

Thank you for purchasing the OptiSmile Water Flosser. This manual will provide instructions for effective and safe use of the product. Please read it carefully before using this device. Store the manual in a safe place for future reference.

Table of contents

1. Package contents
2. Safety instructions
3. Description of parts
4. Operation
5. Maintenance
6. Troubleshooting
7. Disposal

1. Package contents


- 1x Water Flosser
- 1x Charging station
- 2x Nozzles

2. Safety instructions

- For intended use only.
- This is not a toy. Keep out of reach of small children.
- The product is safe to use by children from the age of 8 and by persons with reduced physical or mental capabilities, or with lack of experience or knowledge, but only under supervision and with taking into account any hazards.
- Cleaning and maintenance should be done by or under the supervision of an adult.
- Do not attempt to repair, disassemble or alter the product yourself.
- Do not place, store or use the product near water.
- Do not fill the water tank with water warmer than 40°C or with mouthwash, to avoid damage to the product.
- Do not operate the device without the nozzle.
- Do not hold the nozzle too close to your teeth and gums.
- Stop using the product and contact your doctor, if you experience any discomfort or pain.

3. Description of parts

- Nozzle release button
- ON/OFF button
- Mode button
- Mode indicator
- Handle
- Charging station
- Nozzle
- Water tube
- Filter
- Water tank
- Water tank cap

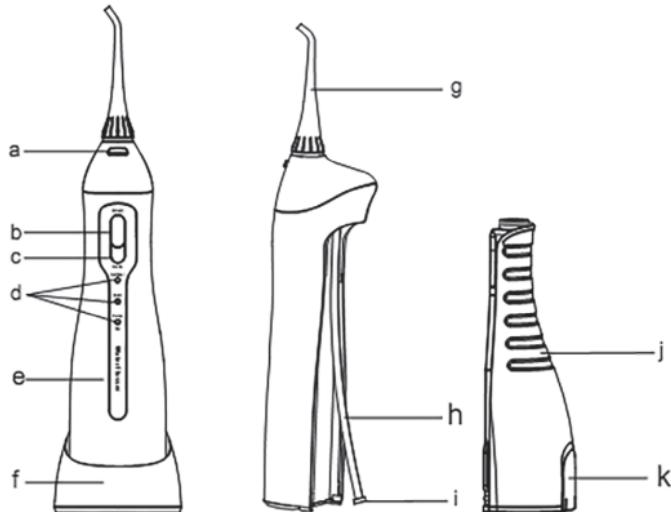Note: we always test the water jet of each product. That is why there may be some distilled water left inside. This is completely safe and will not effect your health.

4. Operation

Charging (Fig. 1)

Make sure that the water flosser is charged before using it for the first time.

1. Plug the adapter of the charging station into the power socket.
2. Place the water flosser on the charging station.
3. The bottom LED light of the mode indicator flashes while charging and stops flashing when the device is fully charged.

Note: it takes 24 hours to fully charge the battery, but the water flosser can be used before that. When the mode indicator flashes during use, the water flosser needs to be charged.

Using the water flosser

1. Place the nozzle on the water flosser (Fig. 2).
2. Open the water tank cap to fill the water flosser with water. Close the cap (Fig. 3). The water tank is detachable for easy refilling (Fig. 4).
3. Press the mode button to select the mode of your choice. You can switch modes while using the water flosser (Fig. 5).
- Normal: for standard cleaning of the mouth
- Soft: for soft cleaning in case of sensitive gums
- Pulsing: for cleaning and massaging the gums
4. Keep the water flosser straight and with the nozzle towards your mouth (90° angle). Guide the water flosser along your gums (Fig. 6).
5. Press the ON/OFF button to operate the water flosser.

After each use (Fig. 8)

1. Open the cap and empty the water tank.
2. Turn the water flosser on to dispose of any remaining water.
3. Turn the device off and close the water tank cap.
4. Wipe the water flosser dry with a clean towel.

Fig. 1

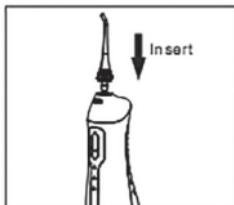
Fig. 2

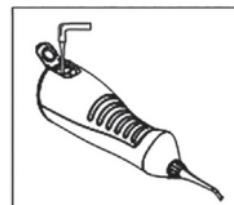
Fig. 3

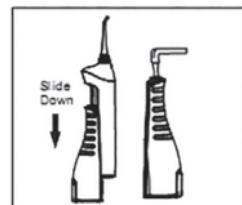
Fig. 4

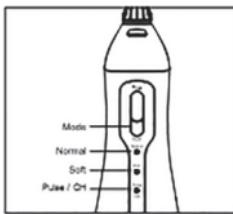
Fig. 5

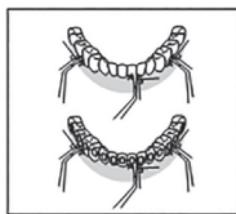
Fig. 6

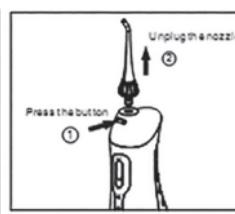
Fig. 7

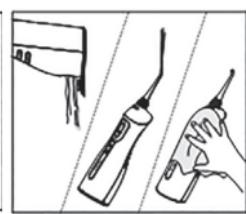
Fig. 8

5. Maintenance

- Do not use water warmer than 40°C.
- Do not immerse the device in water.
- Only use water or a mild non-abrasive detergent to clean the water flosser and all its parts.
- Remove any remaining water (tank, tube or nozzle) to prevent bacteria from spreading.
- After cleaning, wipe dry with a clean and soft cloth.
- When necessary, press the release button to detach the nozzle (Fig. 7). Replace the nozzle every 6 months.
- Store the product in a cool and dry place, and away from direct sunlight.

6. Troubleshooting

Problem	Cause of problem	Solution
The water flosser doesn't turn on.	The product is new or you haven't used it in 3 months.	Charge the water flosser for 24 hours.
The water flosser only works for a couple of minutes.	The battery died.	Please contact our customer service.
	The water flosser isn't charged.	Charge the water flosser for 24 hours.
The water pressure is very low.	The nozzle is congested.	Replace the nozzle.
	The nozzle is bent.	Clean the filter.
	There's dirt in the filter.	
There's no water coming out of the water flosser.	The water tank is empty.	Fill the water tank with water.
	The water flosser is tilted too much.	Keep the water flosser straight.
Water drips from the product and water tank.	The water flosser is tilted too much, causing water to drain out of the vent or the water tank.	Keep the water flosser straight.

7. Disposal

The crossed out wheelie bin symbol on this product is to remind you that electrical products should not be disposed of with household waste. Please dispose of it responsibly at an approved waste or recycling facility. Check with your local authority for recycling advice.

Dispose of all packaging materials in an environmentally friendly manner.

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